Edwards

# Edwards Transfemoral Balloon Catheter

## Catéter balón transfemoral Edwards

## Cateter-balão Transfemoral Edwards

### DIRECTORY

|                                                                       |     |
|-----------------------------------------------------------------------|-----|
| English .....                                                         | 1   |
| Español .....                                                         | 3   |
| Português .....                                                       | 4   |
| Symbol Legend / Inscripción de símbolos /<br>Legenda do Símbolo ..... | 6–7 |

### English

## Instructions for Use

### 1.0 Device Description

The Edwards Transfemoral Balloon Catheter consists of a shaft and balloon with two radiopaque marker bands that indicate the working length of the balloon. The proximal end of the device has a “Y-connector” with a balloon inflation port labeled as “BALLOON” and a guidewire lumen port labeled as “WIRE”.


The inflation parameters are as follows:

**Table 1: Inflation Parameters**

| Model    | Nominal          |                  |                    |
|----------|------------------|------------------|--------------------|
|          | Balloon Diameter | Inflation Volume | Inflation Pressure |
| 9350BC20 | 20 mm            | 16 mL            | 4 atm (405 kPa)    |
| 9350BC23 | 23 mm            | 21 mL            | 4 atm (405 kPa)    |
| 9350BC25 | 25 mm            | 26 mL            | 4 atm (405 kPa)    |

### Edwards Transfemoral Balloon Catheter

THV204

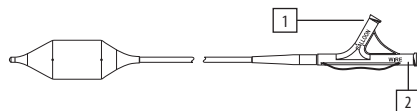Black dots indicate position of radiopaque marker bands.

1 – Balloon Inflation Port

2 – Guidewire Lumen Port

Edwards, Edwards Lifesciences, and the stylized E logo are trademarks of Edwards Lifesciences Corporation. All other trademarks are the property of their respective owners.

Device compatibility specifications are as follows:

**Table 2: Device Compatibility**

| Model    | Max. Guidewire Diameter | Min. Sheath Compatibility |
|----------|-------------------------|---------------------------|
| 9350BC20 | 0.035" (0.89 mm)        | 14F (4.7 mm)              |
| 9350BC23 | 0.035" (0.89 mm)        | 14F (4.7 mm)              |
| 9350BC25 | 0.035" (0.89 mm)        | 16F (5.3 mm)              |

**NOTE:** For proper volume sizing, the balloon catheter should be used with the inflation device provided by Edwards Lifesciences.

### 2.0 Indications

The balloon catheter is indicated for dilation of stenotic native aortic valve leaflets.

### 3.0 Contraindications

The device is contraindicated for patients with:

- evidence of intracardiac mass, thrombus, vegetation, active infection or endocarditis;
- inability to tolerate anticoagulation/antiplatelet therapy.

### 4.0 Warnings

- The device is designed, intended, and distributed for single use only. Do not resterilize or reuse the device. There are no data to support the sterility, nonpyrogenicity, and functionality of the device after reprocessing.
- Observation of the pacing lead throughout the procedure is essential to avoid the potential risk of pacing lead perforation.
- Use only appropriate balloon inflation medium. Do not use air or gaseous medium to inflate the balloon.
- The balloon inflation diameter must not be significantly greater than the annulus diameter being pre-dilated.
- The device is not intended for post-dilatation of deployed transcatheter heart valves.
- While exposed within the body, device advancement and retrieval should not be done without the aid of fluoroscopy. Do not advance or retract the device unless the balloon is fully deflated under vacuum.
- Do not mishandle the balloon catheter or use it if the packaging or any components are not sterile, have been opened or are damaged (i.e. kinked or stretched), or the expiration date has elapsed.

## 5.0 Precautions

The safety and effectiveness of the balloon catheter has not been established in patients who have a congenital unicuspid or congenital bicuspid aortic valve.

## 6.0 Potential Adverse Events

Complications associated with standard catheterization, balloon aortic valvuloplasty, and the use of angiography include, but are not limited to, allergic reaction to anesthesia or to contrast media, injury including perforation or dissection of vessels, thrombosis, emboli formation, renal failure; renal insufficiency, and plaque dislodgement which may result in myocardial infarction, stroke, and/or death. Additional complications may include arrhythmia development, cardiac perforation, conduction system injury, hematoma, infundibulum injury, annular tear or rupture and/or valvular tearing or trauma.

## 7.0 Directions for Use

Dilate native valve leaflets using standard technique and rapid cardiac pacing.

| Step | Procedure                                                                                                                                                                                                                                                                             |
|------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 1    | Prepare vascular access site for catheter insertion and position guidewire using standard techniques.                                                                                                                                                                                 |
| 2    | With the balloon cover in place, flush the guidewire lumen of the Edwards transfemoral balloon catheter with heparinized saline. Attach a high pressure 3-way stopcock to the balloon inflation port.                                                                                 |
| 3    | Prepare a syringe with diluted contrast solution (15:85 medium to saline dilution) and attach to the stopcock.                                                                                                                                                                        |
| 4    | Fill the inflation device with excess diluted contrast medium relative to the indicated volume, attach in the locked position to the stopcock, close the stopcock to the inflation device.                                                                                            |
| 5    | Slowly pull vacuum with the syringe repeatedly to remove air, leaving zero pressure in the system.                                                                                                                                                                                    |
| 6    | Close the stopcock to the balloon catheter. Gradually remove contrast medium into the syringe to achieve the appropriate volume (as specified in Table 1: Inflation Parameters). Lock the inflation device, close the stopcock to the syringe and remove the syringe from the system. |
| 7    | Remove the balloon cover and hydrate the length of the balloon catheter.                                                                                                                                                                                                              |
| 8    | Advance the balloon catheter over the guidewire, through the introducer sheath, across the aortic valve, and position the balloon markers at the intended site.                                                                                                                       |
| 9    | Ensure hemodynamic stability is established and begin rapid pacing. Once the blood pressure has decreased to 50 mmHg or below, balloon inflation can commence.                                                                                                                        |
| 10   | Fully and rapidly inflate the balloon with the inflation device. In case of balloon instability, repeat balloon inflation while ensuring rapid ventricular pacing. When the balloon is fully deflated, the pacing should be turned off.                                               |

## 8.0 How Supplied

Supplied pouched and sterilized by ethylene oxide.

## 9.0 Storage

Store in a cool, dry place.

## 10.0 Device Disposal

Used devices may be handled and disposed of in the same manner as hospital waste and biohazardous materials. There are no special risks related to the disposal of these devices.

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© Copyright 2018, Edwards Lifesciences LLC
All rights reserved
**Edwards Lifesciences Services GmbH**
**Edwards Lifesciences LLC**