Edwards

Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve System

Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve

Edwards COMMANDER Delivery System

Edwards SAPIEN 3 og SAPIEN 3 Ultra transkateter hjerteklaffsystem

Edwards SAPIEN 3 og SAPIEN 3 Ultra transkateter hjerteklaff

Edwards COMMANDER leveringssystem

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English

Instructions for Use

Implantation of the transcatheter heart valves should be performed only by physicians who have received Edwards Lifesciences training. The implanting physician should be experienced in standard catheterization techniques. It is at the physician's discretion to choose the appropriate access route to implant the THV based on the patient anatomy and associated risks.

1.0 Device Description

Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve System

The Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve (THV) system consists of the Edwards SAPIEN 3 and SAPIEN 3 Ultra transcatheter heart valves and delivery systems.

• Edwards SAPIEN 3 Transcatheter Heart Valve (Figure 1)

The Edwards SAPIEN 3 transcatheter heart valve (THV) is comprised of a balloon-expandable, radiopaque, cobalt-chromium frame, trileaflet bovine pericardial tissue valve, and polyethylene terephthalate (PET) inner and outer fabric skirts. The leaflets are treated according to the Carpentier-Edwards ThermaFix process.

• Edwards SAPIEN 3 Ultra Transcatheter Heart Valve (Figure 2)

The Edwards SAPIEN 3 Ultra transcatheter heart valve (THV) is comprised of a balloon-expandable, radiopaque, cobalt-chromium frame, trileaflet bovine pericardial tissue valve, and polyethylene terephthalate (PET) inner and outer fabric skirts. The leaflets are treated according to the Carpentier-Edwards ThermaFix process.

Edwards, Edwards Lifesciences, the stylized E logo, Carpentier-Edwards, Edwards COMMANDER, Edwards eSheath, Edwards SAPIEN, Edwards SAPIEN 3, Edwards SAPIEN 3 Ultra, eSheath, INSPIRIS, INSPIRIS RESILIA, Qualcrimp, RESILIA, SAPIEN, SAPIEN 3, SAPIEN 3 Ultra, ThermaFix, and VFit are trademarks of Edwards Lifesciences Corporation. All other trademarks are the property of their respective owners.

The THV is intended to be implanted in a native annulus size range associated with the three-dimensional area of the aortic annulus measured at the basal ring during systole as provided in Table 1:

Table 1

Native Valve Annulus Size (TEE)*	Native Valve Annulus Size (CT)		THV Size
	Area	Area Derived Diameter	
16 – 19 mm	273 – 345 mm ²	18.6 – 21.0 mm	20 mm
18 – 22 mm	338 – 430 mm ²	20.7 – 23.4 mm	23 mm
21 – 25 mm	430 – 546 mm ²	23.4 – 26.4 mm	26 mm
24 – 28 mm	540 – 683 mm ²	26.2 – 29.5 mm	29 mm

THV size recommendations are based on native valve annulus size, as measured by transesophageal echocardiography (TEE) or computed tomography (CT). Patient anatomical factors and multiple imaging modalities should be considered during THV size selection.

Note: Risks associated with undersizing and oversizing should be considered to minimize the risk of paravalvular leak, migration, and/or annular rupture.

*Due to limitations in two-dimensional images, 2-D TEE imaging should be supplemented with 3-D area measurements.

Sizing recommendations for implanting the THV in a failing bioprosthesis are provided in Table 2:

Table 2

Surgical Valve True Inner Diameter (ID) ^[1]	THV-in-THV (Native Valve Annulus Size)	THV Size
16.5 – 19.0 mm	18.6 – 21.0 mm	20 mm
18.5 – 22.0 mm	20.7 – 23.4 mm	23 mm
22.0 – 25.0 mm	23.4 – 26.4 mm	26 mm
25.0 – 28.5 mm	26.2 – 29.5 mm	29 mm

Note: Surgical valve 'True ID' may be smaller than the labeled valve size. For THV-in-THV, the native valve annulus size should be considered to determine the appropriate THV size to implant. For a failing stentless bioprosthesis, consider sizing recommendations for a native annulus. The dimensions of the failed bioprosthesis should be determined so that the appropriate THV size can be implanted; and is best determined by using computed tomography, magnetic resonance imaging, and/or transesophageal echocardiography.

Sizing recommendations for implanting the THV in a failing INSPIRIS RESILIA aortic surgical bioprosthesis in sizes 19 – 25 mm, based on bench testing, are provided in Table 3 below:

Table 3

INSPIRIS RESILIA Aortic Valve (model 11500A)* Labeled Size	THV Size
19 mm	20 mm or 23 mm
21 mm	23 mm or 26 mm
23 mm	26 mm
25 mm	29 mm

*INSPIRIS RESILIA aortic valve model 11500A sizes 19 – 25 mm incorporate VFit technology which consists of expandable bands and fluoroscopically visible size markers designed for potential future valve-in-valve procedures. Clinical data are not currently available on the INSPIRIS RESILIA aortic valve Model 11500A valve-in-valve procedure or expansion feature. The impact of tissue ingrowth on the expansion feature of the INSPIRIS RESILIA aortic valve has not been assessed.

WARNING: Do not perform stand-alone balloon aortic valvuloplasty procedures in the INSPIRIS RESILIA aortic valve for the sizes 19 – 25 mm. This may expand the valve causing aortic incompetence, coronary embolism or annular rupture.

Note: INSPIRIS RESILIA aortic valve model 11500A sizes 27 – 29 mm do not incorporate VFit technology and therefore follow the surgical valve True ID sizing provided in Table 2.

Note: Exact volume required to deploy the THV may vary depending on the bioprosthesis inner diameter. Factors such as calcification and pannus tissue growth may not be accurately visualized in imaging and may reduce the effective inner diameter of the failing bioprosthesis to a size smaller than the 'True ID'.

These factors should be considered and assessed in order to determine the most appropriate THV size to achieve nominal THV deployment and sufficient anchoring. Do not exceed the rated burst pressure. See inflation parameters in Table 4.

• Edwards COMMANDER Delivery System (Figure 4)

The Edwards COMMANDER delivery system facilitates the placement of the bioprosthesis.

It consists of a Flex Catheter to aid in valve alignment to the balloon, tracking, and positioning of the THV. The delivery system includes a tapered tip to facilitate crossing of the valve. The handle contains a Flex Wheel to control flexing of the Flex Catheter, and a Balloon Lock and Fine Adjustment Wheel to facilitate valve alignment and positioning of the valve within the target location. A stylet is included within the guidewire lumen of the delivery system. The Balloon Catheter has radiopaque Valve Alignment Markers defining the working length of the balloon. A radiopaque Center Marker in the balloon is provided to help with valve positioning. A radiopaque Triple Marker proximal to the balloon indicates the Flex Catheter position during deployment.

The inflation parameters for valve deployment are:

Table 4

Model	Nominal Balloon Diameter	Nominal Inflation Volume	Rated Burst Pressure (RBP)
9610TF20	20 mm	11 ml	7 atm
9610TF23	23 mm	17 ml	7 atm
9610TF26	26 mm	23 ml	7 atm
9610TF29	29 mm	33 ml	7 atm

• Edwards Sheath

Refer to the sheath instructions for use for device description.

• Qualcrimp Crimping Accessory

The Qualcrimp crimping accessory is used during THV crimping (Figure 3).

• Loader

The loader is used to aid insertion of the delivery system into the sheath (Figure 5).

• Edwards Crimper and Crimp Stopper (Figure 6)

The Edwards crimper reduces the diameter of the valve to mount it onto the delivery system. The crimper is comprised of a housing and a compression mechanism that is closed with a handle located on the housing. A 2-piece crimp stopper is used to crimp the valve to its intended diameter.

• Inflation Device

An inflation device with locking mechanism is used during valve deployment.

Note: For proper volume sizing, the delivery system must be used with the inflation device provided by Edwards Lifesciences.

2.0 Intended Use

The bioprosthesis is intended for use in patients requiring heart valve replacement. The delivery system and accessories are intended to facilitate the placement of the bioprosthesis via the transfemoral, transseptal, subclavian/axillary access approaches.

3.0 Indications

- The Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve system is indicated for use in patients with heart disease due to native calcific aortic stenosis at any or all levels of surgical risk for open heart surgery.
- The Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve system is indicated for patients with symptomatic heart disease due to failure (stenosed, insufficient or combined) of an aortic transcatheter bioprosthetic or surgical aortic or mitral bioprosthetic valve who are judged by a heart team, including a cardiac surgeon, to be at high or greater risk for open surgical therapy (i.e., predicted risk of surgical mortality $\geq 8\%$ at 30 days, based on the Society of Thoracic Surgeons (STS) risk score and other clinical co-morbidities unmeasured by the STS risk calculator).

4.0 Contraindications

Use of the Edwards SAPIEN 3 and SAPIEN 3 Ultra THV system is contraindicated in patients who:

- Cannot tolerate anticoagulation/antiplatelet regimen or who have active bacterial endocarditis or other active infections.

5.0 Warnings

- The devices are designed, intended, and distributed STERILE for single use only. **Do not resterilize or reuse the devices.** There are no data to support the sterility, nonpyrogenicity, and functionality of the devices after reprocessing.
- Correct sizing of the THV is essential to minimize the risk of paravalvular leak, migration, and/or annular rupture.
- The physician must verify correct orientation of the THV prior to its implantation.
- Accelerated deterioration of the THV may occur in patients with an altered calcium metabolism.
- Observation of the pacing lead throughout the procedure is essential to avoid the potential risk of pacing lead perforation.
- The THV must remain hydrated at all times and cannot be exposed to solutions, antibiotics, chemicals, etc. other than its shipping storage solution and sterile physiologic saline solution to prevent leaflet damage that may impact valve functionality. THV leaflets mishandled or damaged during any part of the procedure will require replacement of the THV.
- Patients with hypersensitivities to cobalt, nickel, chromium, molybdenum, titanium, manganese, silicon, bovine tissue, and/or polymeric materials may have an allergic reaction to these materials.
- Do not use the THV if the tamper evident seal is broken, as sterility may be compromised.
- Do not use the THV if the temperature indicator has been activated, as valve function may be compromised.
- Do not use the THV if the expiration date has elapsed, as either sterility or valve function may be compromised.
- Do not use the THV if the storage solution does not completely cover the THV or the THV is damaged.
- Do not mishandle the delivery system or use the delivery system and accessory devices if the packaging sterile barriers and any components have been opened or damaged, cannot be flushed, or the expiration date has elapsed.

- Patient injury could occur if the delivery system is not un-flexed prior to removal.
- Access characteristics such as severe obstructive or circumferential calcification, severe tortuosity, vessel diameters less than 5.5 mm (for size 20, 23 and 26 mm SAPIEN 3/SAPIEN 3 Ultra transcatheter heart valve) or 6.0 mm (for 29 mm SAPIEN 3 transcatheter heart valve) may preclude safe placement of the sheath and should be carefully assessed prior to the procedure.
- Valve recipients should be maintained on anticoagulant/antiplatelet therapy, except when contraindicated, to minimize the risk of valve thrombosis or thromboembolic events, as determined by their physicians. This device has not been tested for use without anticoagulation.
- The procedure should be conducted under fluoroscopic guidance. Some fluoroscopically guided procedures are associated with a risk of radiation injury to the skin. These injuries may be painful, disfiguring, and long-lasting.
- Caution should be exercised in implanting a valve in patients with clinically significant coronary artery disease.
- Patients with pre-existing bioprostheses should be carefully assessed prior to implantation of the valve to ensure proper valve positioning and deployment.
- Balloon valvuloplasty should be avoided in the treatment of failing bioprostheses as this may result in embolization of bioprosthesis material and mechanical disruption of the valve leaflets.

6.0 Precautions

- Long-term durability has not been established for the THV. Regular medical follow-up is advised to evaluate valve performance.
- Glutaraldehyde may cause irritation of the skin, eyes, nose and throat. Avoid prolonged or repeated exposure to, or breathing of, the solution. Use only with adequate ventilation. If skin contact occurs, immediately flush the affected area with water; in the event of contact with eyes, seek immediate medical attention. For more information about glutaraldehyde exposure, refer to the Material Safety Data Sheet available from Edwards Lifesciences.
- The safety and effectiveness of the THV implantation has not been established in patients who have:
 - Non-calcified aortic annulus
 - Congenital unicuspid aortic valve
 - Pre-existing prosthetic ring in any position
 - Severe ventricular dysfunction with ejection fraction < 20%
 - Hypertrophic cardiomyopathy with or without obstruction
 - Aortic stenosis characterized by a combination of AV low flow, low gradient
 - A concomitant paravalvular leak where the failing bioprosthesis is not securely fixed in the native annulus or is not structurally intact (e.g., wireframe fracture)
 - A partially detached leaflet of the failing bioprosthesis that in the aortic position may obstruct a coronary ostium
 - Bulky calcified aortic valve leaflets in close proximity to coronary ostia.
 - Significant aortic disease, including abdominal aortic or thoracic aneurysm defined as maximal luminal diameter 5 cm or greater; marked tortuosity (hyperacute bend), aortic arch atheroma (especially if thick (> 5 mm), protruding, or ulcerated) or narrowing (especially with calcification and surface irregularities) of the abdominal or thoracic aorta, severe "unfolding" and tortuosity of the thoracic aorta.
- If a significant increase in resistance occurs when advancing the catheter through the vasculature, stop advancement and investigate the cause of resistance before proceeding. Do not force passage, as this could increase the risk of vascular complications. As compared to SAPIEN 3, system advancement force may be higher with the use of SAPIEN 3 Ultra THV in tortuous/challenging vessel anatomies.
- Appropriate antibiotic prophylaxis is recommended post-procedure in patients at risk for prosthetic valve infection and endocarditis.
- Additional precautions for transseptal replacement of a failed mitral valve bioprosthesis include: presence of devices or thrombus or other abnormalities in the caval vein precluding safe transvenous femoral access for transseptal approach; presence of Atrial Septal Occluder Device or calcium preventing safe transseptal access.
- Special care must be exercised in mitral valve replacement if chordal preservation techniques were used in the primary implantation to avoid entrapment of the subvalvular apparatus.

- Based on the treating physician's consideration of risks and benefits, the valve may be implanted in relatively young patients, although the longer-term durability is still the subject of ongoing clinical research.
- Do not overinflate the deployment balloon, as this may prevent proper valve leaflet coaptation and thus impact valve functionality.
- Patients with pre-existing mitral valve devices should be carefully evaluated before implantation of the THV to ensure proper THV positioning and deployment.
- The risks of subclavian/axillary access are low and acceptable, but subclavian/axillary access should be considered when the physician determines there is an increased risk associated with transfemoral access.
- For left axillary approach, a left subclavian takeoff angle $\sim \geq 90^\circ$ from the aortic arch causes sharp angles, which may be responsible for potential sheath kinking, subclavian/axillary dissection and aortic arch damage.
- For left/right axillary approach, ensure there is flow in the Left Internal Mammary Artery (LIMA) / Right Internal Mammary Artery (RIMA) respectively during procedure and monitor pressure in homolateral radial artery.
- Residual mean gradient may be higher in a "THV-in-failing bioprosthesis" configuration than that observed following implantation of the valve inside a native aortic annulus using the same size device. Patients with elevated mean gradient post procedure should be carefully followed. It is important that the manufacturer, model and size of the preexisting bioprosthetic valve be determined, so that the appropriate valve can be implanted and a prosthesis-patient mismatch be avoided. Additionally, pre-procedure imaging modalities must be employed to make as accurate a determination of the inner diameter as possible.

7.0 Potential Adverse Events

Potential risks associated with the overall procedure including access, cardiac catheterization, local and/or general anesthesia:

- Death
- Stroke/transient ischemic attack, clusters or neurological deficit
- Paralysis
- Permanent disability
- Respiratory insufficiency or respiratory failure
- Hemorrhage requiring transfusion or intervention
- Cardiovascular injury including perforation or dissection of vessels, ventricle, atrium, septum, myocardium or valvular structures that may require intervention
- Pericardial effusion or cardiac tamponade
- Thoracic bleeding
- Embolization including air, calcific valve material or thrombus
- Infection including septicemia and endocarditis
- Heart failure
- Myocardial ischemia or infarction
- Renal insufficiency or renal failure
- Conduction system defect which may require a permanent pacemaker
- Arrhythmias including ventricular fibrillation (VF) and ventricular tachycardia (VT)
- Retroperitoneal bleed
- Arteriovenous (AV) fistula or pseudoaneurysm
- Reoperation
- Ischemia or nerve injury or brachial plexus injury or compartment syndrome
- Restenosis
- Pulmonary edema
- Pleural effusion
- Bleeding, bleeding requiring transfusion or intervention
- Anemia
- Vessel thrombosis/occlusion
- Abnormal lab values (including electrolyte imbalance)
- Hypertension or hypotension
- Allergic reaction to anesthesia, contrast media, device materials or bovine pericardial tissue
- Hematoma
- Syncope
- Pain or changes (e.g., wound infection, hematoma, and other wound care complications) at the access site
- Exercise intolerance or weakness

- Inflammation
- Angina
- Vasovagal response
- Heart murmur
- Fever

Additional potential risks associated with the TAVR procedure, the bioprosthesis, and the use of its associated devices and accessories include:

- Cardiac arrest
- Cardiogenic shock
- Emergency cardiac surgery
- Cardiac failure or low cardiac output
- Coronary flow obstruction/transvalvular flow disturbance
- Device thrombosis requiring intervention
- Valve thrombosis
- Device embolization
- Device migration or malposition requiring intervention
- Left ventricular outflow tract obstruction
- Valve deployment in unintended location
- Valve stenosis
- Vessel spasm
- Structural valve deterioration (wear, fracture, calcification, leaflet tear/tearing from the stent posts, leaflet retraction, suture line disruption of components of the prosthetic valve, thickening, stenosis)
- Device degeneration
- Paravalvular or transvalvular leak
- Valve regurgitation
- Hemolysis
- Device explants
- Nonstructural dysfunction
- Mechanical failure of delivery system, and/or accessories, including balloon rupture and tip separation
- Non-emergent reoperation
- Allergic/immunologic reaction to the implant
- Injury to mitral valve

For a patient/user/third party in the European Economic area; if, during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and your national competent authority, which can be found at https://ec.europa.eu/growth/sectors/medical-devices/contacts_en.

8.0 Directions for Use

8.1 System Compatibility

Table 5

Product Name	20 mm System	23 mm System	26 mm System	29 mm System
	Model			
Edwards SAPIEN 3 Transcatheter Heart Valve	9600TFX (20 mm)	9600TFX (23 mm)	9600TFX (26 mm)	9600TFX (29 mm)
Edwards COMMANDER Delivery System	9610TF20	9610TF23	9610TF26	9610TF29
Edwards eSheath Introducer Set	9610ES14			9610ES16
Inflation Device	96402			96406
Edwards Crimper	9600CR			
Qualcrimp Crimping Accessory, Crimp Stopper and Loader provided by Edwards Lifesciences				

Table 6

Product Name	20 mm System	23 mm System	26 mm System
	Model		
Edwards SAPIEN 3 Ultra Transcatheter Heart Valve	9750TFX (20 mm)	9750TFX (23 mm)	9750TFX (26 mm)
Edwards COMMANDER Delivery System	9610TF20	9610TF23	9610TF26
Edwards eSheath Introducer Set	9610ES14		
Inflation Device	96402		
Edwards Crimper	9600CR		
Qualcrimp Crimping Accessory, Crimp Stopper and Loader provided by Edwards Lifesciences			

Additional Equipment:

- Balloon catheter per the discretion of the physician
- 20 cc syringe or larger
- 50 cc syringe or larger
- High-pressure 3-way stopcock
- Standard cardiac catheterization lab equipment and supplies, and access to standard heart valve operating room equipment and supplies
- Fluoroscopy (fixed, mobile or semi-mobile fluoroscopy systems appropriate for use in percutaneous coronary interventions)
- Transesophageal or transthoracic echocardiography capabilities
- Exchange length 0.035 in (0.89 mm) extra-stiff guidewire
- Temporary pacemaker (PM) and pacing lead
- Sterile rinsing basins, physiological saline, heparinized saline, 15% diluted radiopaque contrast medium
- Sterile table for THV and accessories preparation

8.2 THV Handling and Preparation

Follow sterile technique during device preparation and implantation.

8.2.1 THV Rinsing Procedure

Before opening the valve jar, carefully examine for evidence of damage (e.g., a cracked jar or lid, leakage, or broken or missing seals).

CAUTION: If the container is found to be damaged, leaking, without adequate sterilant, or missing intact seals, the THV must not be used for implantation, as sterility may be compromised.

1. Set up two (2) sterile bowls with at least 500 ml of sterile physiological saline to thoroughly rinse the THV.
2. Carefully remove the valve/holder assembly from the jar without touching the tissue. Verify the valve serial identification number with the number on the jar lid and record in the patient information documents. Inspect the valve for any signs of damage to the frame or tissue.
3. Rinse the THV as follows:
 - a) Place the THV in the first bowl of sterile, physiological saline. Be sure the saline solution completely covers the THV and holder.
 - b) With the valve and holder submerged, slowly agitate (to gently swirl the valve and holder) back and forth for a minimum of 1 minute.
 - c) Transfer the THV and holder to the second rinsing bowl of sterile physiological saline and gently agitate for at least one more minute. Ensure the rinse solution in the first bowl is not used.
 - d) The valve should be left in the final rinse solution until needed to prevent the tissue from drying.

CAUTION: Do not allow the valve to come into contact with the bottom or sides of the rinse bowl during agitation or swirling in the rinse solution. Direct contact between the identification tag and valve is also to be avoided during the rinse procedure. No other objects should be placed in the rinse bowls. The valve should be kept hydrated to prevent the tissue from drying.

8.2.2 Prepare the System

1. Visually inspect all the components for damage. Ensure the delivery system is fully unflexed and the balloon catheter is fully advanced in the flex catheter.

WARNING: To prevent possible damage to the balloon shaft, ensure that the proximal end of the balloon shaft is not subjected to bending.

2. Flush the delivery system with heparinized saline through the flush port.
3. Remove the distal balloon cover from the delivery system. Remove the stylet from the distal end of the guidewire lumen and set aside.
4. Flush the guidewire lumen with heparinized saline. Insert the stylet back into the guidewire lumen.

Note: Failure to replace the stylet in the guidewire lumen may result in damage to the lumen during the THV crimping process.

5. Place the delivery system into the Default Position (end of strain relief is aligned between the two white markers on the balloon shaft) and make sure that the flex catheter tip is covered by the proximal balloon cover.
6. Unscrew the loader cap from the loader and flush the loader cap with heparinized saline.
7. Place the loader cap onto the delivery system with the inside of the cap oriented towards the distal tip.
Fully advance the balloon catheter in the flex catheter.
Peel off the proximal balloon cover over the blue section of the balloon shaft.
8. Attach a 3-way stopcock to the balloon inflation port. Fill a 50 cc or larger syringe with 15-20 ml of diluted contrast medium and attach to the 3-way stopcock.
9. Fill the inflation device with excess volume of diluted contrast medium relative to the indicated inflation volume. Lock and attach to the 3-way stopcock. Close stopcock to the inflation device.
10. Pull vacuum with the syringe to remove air. Slowly release the plunger to ensure that the contrast medium enters the lumen of the delivery system. Repeat until all air bubbles are removed from the system. Leave zero-pressure in the system.

WARNING: Ensure there is no residual fluid left in the balloon to avoid potential difficulty with valve alignment during the procedure.

Close stopcock to the delivery system.

11. Rotate the knob of the inflation device to remove the contrast medium into the syringe and achieve the appropriate volume required to deploy the THV.
Close the stopcock to the syringe and remove syringe.
12. Verify that the inflation volume in the inflation device is correct.

CAUTION: Maintain the inflation device in the locked position until THV deployment to minimize the risk of premature balloon inflation and subsequent improper THV deployment.

8.2.3 Mount and Crimp the THV onto the Delivery System


1. Completely submerge the Qualcrimp crimping accessory in a bowl of 100 ml physiological saline. Gently compress until fully saturated. Swirl for a minimum of 1 minute. Repeat this process in a second bowl.
2. Remove the THV from the holder and remove the ID tag.
3. Rotate the crimper handle until the aperture is fully open. Attach the 2-piece Crimp Stopper to the base of the crimper and click into place.
4. If necessary, partially crimp the THV in the crimper until it fits snugly inside the Qualcrimp crimping accessory.

Note: Partial crimping is not necessary for the 20 mm valve.

5. Place the Qualcrimp crimping accessory over the THV aligning the edge of the Qualcrimp crimping accessory with the outflow of the THV.
6. Place the THV and Qualcrimp crimping accessory in crimper aperture. Insert the delivery system coaxially within the THV 2-3 mm distal to the blue balloon shaft (in the Valve Crimp Section) of the delivery system with the orientation of the valve on the delivery system as described below:


Antegrade approach:

Inflow (outer skirt end) of the THV towards the **proximal end** of the delivery system.

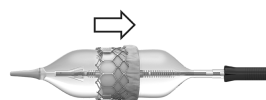Retrograde approach:

Inflow (outer skirt end) of the THV towards the **distal end** of the delivery system.

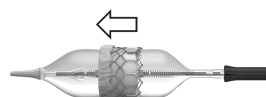7. Center the balloon shaft coaxially within the THV. Crimp the THV until it reaches the Qualcrimp stop.
8. Remove the Qualcrimp crimping accessory from the THV and Qualcrimp stop from the Crimp Stopper, leaving the Final Stop in place.
9. Center the THV within the crimper aperture. Fully crimp the THV until it reaches the Final Stop and hold for 5 seconds. Repeat this crimp step two (2) more times for a total of 3 crimps.

Note: Ensure that the Valve Crimp Section is coaxial within the THV.

10. Pull the balloon shaft and engage the Balloon Lock so the delivery system is in Default Position.
11. Flush the loader with heparinized saline. Immediately advance the THV into the loader until the tapered tip of the delivery system is exposed.

CAUTION: The THV should not remain fully crimped and/or in the loader for over 15 minutes, as leaflet damage may result and impact valve functionality.

12. Attach the loader cap to the loader, re-flush the Flex Catheter and close the stopcock to the delivery system. Remove the stylet and flush the guidewire lumen of the delivery system.

CAUTION: Keep THV hydrated until ready for implantation to prevent damage to the leaflets which may impact valve functionality.

WARNING: The physician must verify correct orientation of the THV prior to its implantation to prevent the risk of severe patient harm.

8.3 Native Valve Predilation and THV Delivery

Native valve predilation and THV delivery should be performed under local and/or general anesthesia with hemodynamic monitoring in a catheterization lab/hybrid operating room with fluoroscopic and echocardiographic imaging capabilities.

Administer heparin to maintain the ACT at ≥ 250 sec.

CAUTION: Contrast media usage should be monitored to reduce the risk of renal injury.

Balloon valvuloplasty should be avoided in the treatment of failing bioprostheses as this may result in embolization of bioprosthesis material and mechanical disruption of the valve leaflets.

8.3.1 Baseline Parameters

1. Perform an angiogram with the projection of the valve perpendicular to the view.
2. For aortic implantation, evaluate the distance of the left and right coronary ostia from the aortic annulus in relation to the THV frame height.
3. Introduce a pacemaker (PM) lead and position appropriately.
4. Set the stimulation parameters to obtain 1:1 capture, and test pacing.

8.3.2 Native Valve Predilation

Pre-dilate the native aortic valve, per the discretion of the physician, according to the instructions for use for the selected balloon aortic valvuloplasty catheter.

CAUTION: Valve implantation should not be carried out if balloon cannot be fully inflated during valvuloplasty.

8.3.3 THV Delivery

1. Prepare the Edwards sheath introducer set per its instructions for use.
2. If necessary, predilate the vessel.
3. Introduce the sheath per its instructions for use.
4. Insert the loader assembly into the sheath until the loader stops.
5. Advance the delivery system until the THV exits the sheath.

CAUTION: For iliofemoral access, the THV should not be advanced through the sheath if the sheath tip is not past the bifurcation to minimize the risk of vessel damage.

CAUTION: The THV should not remain in the sheath for over 5 minutes as leaflet damage may result and impact valve functionality.

6. In a straight section of the vasculature, initiate valve alignment by disengaging the Balloon Lock and pulling the balloon catheter straight back until part of the Warning Marker is visible. Do not pull past the Warning Marker.

WARNING: To prevent possible damage to the balloon shaft, ensure that the proximal end of the balloon shaft is not subjected to bending.

WARNING: If valve alignment is not performed in a straight section, there may be difficulties performing this step which may lead to delivery system damage and inability to inflate the balloon. Utilizing alternate fluoroscopic views may help with assessing curvature of the anatomy. If excessive tension is experienced during valve alignment, repositioning the delivery system to a different straight section of the aorta and relieving compression (or tension) in the system will be necessary.

Engage the Balloon Lock.

Utilize the Fine Adjustment Wheel to position the THV between the Valve Alignment Markers.

Note: Do not turn the Fine Adjustment Wheel if the Balloon Lock is not engaged.

WARNING: Do not position the THV past the distal Valve Alignment Marker to minimize the risk of improper THV deployment or THV embolization.

CAUTION: Maintain guidewire position during valve alignment to prevent loss of guidewire position.

7. Utilize the Flex wheel to access and cross the valve.

Note: Verify the orientation of the Edwards logo to ensure proper articulation.

Note: The delivery system articulates in a direction opposite from the flush port.

8. Disengage the Balloon Lock and retract the tip of the Flex Catheter to the center of the Triple Marker. Engage the Balloon Lock.
9. Position the THV with respect to the valve.
10. As necessary, utilize the Flex wheel to adjust the co-axiality of the THV and the Fine Adjustment Wheel to adjust the position of the THV.
11. Before deployment, ensure that the THV is correctly positioned between the Valve Alignment Markers and the Flex Catheter tip is over the Triple Marker.
12. Begin THV deployment:
 - a) Unlock the inflation device provided by Edwards Lifesciences.
 - b) Ensure hemodynamic stability is established and begin rapid pacing; once arterial blood pressure has decreased to 50 mmHg or below, balloon inflation can commence.
 - c) Using slow controlled inflation, deploy the THV with the entire volume in the inflation device, hold for 3 seconds and confirm that the barrel of the inflation device is empty to ensure complete inflation of the balloon.
 - d) Deflate the balloon. When the balloon has been completely deflated turn off the pacemaker.

8.3.4 System Removal

1. Unflex the delivery system while retracting the device. Verify that the Flex Catheter tip is locked over the Triple Marker. Retract the loader to the proximal end of the delivery system. Remove the delivery system from the sheath.

Note: For subclavian-axillary approach, keep delivery system inside sheath until ready to remove all devices as one unit.

CAUTION: Completely unflex the delivery system prior to removal to minimize the risk of vascular injury.

8.4 Verification of THV Position and Measurements

Measure and record hemodynamic parameters.

1. Perform an angiogram to evaluate device performance and coronary patency, where applicable.
2. Measure and record the transvalvular pressure gradients.
3. Remove all devices when the ACT level is appropriate (e.g., reaches < 150 sec).
Refer to the introducer sheath instructions for use for device removal.
4. Close the access site.

9.0 How Supplied

STERILE: The valve is supplied sterilized with glutaraldehyde solution.

The delivery system and accessories are supplied sterilized by ethylene oxide gas.

The THV is supplied nonpyrogenic packaged in buffered glutaraldehyde, in a plastic jar to which a tamper evident seal has been applied. Each jar is shipped in a shelf box containing a temperature indicator to detect exposure of the THV to extreme temperature. The shelf box is enclosed in Styrofoam prior to shipping.

9.1 Storage

The THV must be stored at 10 °C to 25 °C (50 °F to 77 °F). Each jar is shipped in an enclosure containing a temperature indicator to detect exposure of the THV to extreme temperature.

The delivery system and accessories should be stored in a cool, dry place.

10.0 MR Safety

MR Conditional

Non-clinical testing has demonstrated that the Edwards SAPIEN 3 Ultra transcatheter heart valve and Edwards SAPIEN 3 transcatheter heart valve are MR Conditional. A patient with this device can be scanned safely, immediately after placement of this device under the following conditions:

- Static magnetic field of 1.5 tesla (T) or 3.0 tesla (T)
- Maximum spatial gradient field of 2500 Gauss/cm (25 T/m) or less
- Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2.0 W/kg (Normal Operating Mode)

Under the scan conditions defined above, the transcatheter heart valve is expected to produce a maximum temperature rise of 3.0 °C after 15 minutes of continuous scanning.

In non-clinical testing, the image artifact caused by the device extends as far as 14.5 mm from the implant for spin echo images and 30 mm for gradient echo images when scanned in a 3.0 T MRI system. The artifact obscures the device lumen in gradient echo images.

The implant has not been evaluated in MR systems other than 1.5 T or 3.0 T.

For valve-in-valve implantation or in the presence of other implants, please refer to the MRI safety information for the surgical valve or other devices prior to MR imaging.

11.0 Qualitative and Quantitative Information related to the THV

This device contains the following substance(s) defined as CMR 1B in a concentration above 0.1% weight by weight:

Cobalt; CAS No. 7440-48-4; EC No. 231-158-0

Current scientific evidence supports that medical devices manufactured from cobalt alloys or stainless steel alloys containing cobalt do not cause an increased risk of cancer or adverse reproductive effects.

For THV, the following tables show the qualitative and quantitative information on the materials and substances:

Table 7: SAPIEN 3 Transcatheter Heart Valve

Substance	CAS	Model Mass Range (mg)
Cobalt	7440-48-4	131 - 427
Nickel	7440-02-0	148 - 405
Chromium	7440-47-3	85.2 - 230
Polyethylene terephthalate	25038-59-9	102 - 170
Collagens, bovine, polymers with glutaraldehyde	2370819-60-4	58.3 - 141
Molybdenum	7439-98-7	40.3 - 115
Polytetrafluoroethylene	9002-84-0	17.5 - 25.5
Polyethylene	9002-88-4	14.2 - 19.7
Iron	7439-89-6	0 - 10.9
Titanium	7440-32-6	0 - 10.9
Manganese	7439-96-5	0 - 1.64
Silicon	7440-21-3	0 - 1.64
Titanium dioxide	13463-67-7	0.219 - 0.752
Polybutylate	24936-97-8	0.273 - 0.383
Carbon	7440-44-0	0 - 0.274
Antimony trioxide	1309-64-4	0.112 - 0.190
Boron	7440-42-8	0 - 0.164
Phosphorus	7723-14-0	0 - 0.164
Sulfur	7704-34-9	0 - 0.109
D&C Green No. 6	128-80-3	0.0394 - 0.0578
Silicon dioxide	7631-86-9	0.00422 - 0.00592
Erucamide	112-84-5	0.000683 - 0.00128
4-Dodecylbenzenesulfonic acid	121-65-3	0.000286 - 0.000430

Table 8: SAPIEN 3 Ultra Transcatheter Heart Valve

Substance	CAS	Model Mass Range (mg)
Cobalt	7440-48-4	131 - 314
Nickel	7440-02-0	148 - 298
Polyethylene terephthalate	25038-59-9	142 - 212
Chromium	7440-47-3	85.2 - 169
Collagens, bovine, polymers with glutaraldehyde	2370819-60-4	58.3 - 97.5
Molybdenum	7439-98-7	40.3 - 84.6
Polyethylene	9002-88-4	19.4 - 22.0
Polytetrafluoroethylene	9002-84-0	12.3 - 15.1
Iron	7439-89-6	0 - 8.06
Titanium	7440-32-6	0 - 8.06
Manganese	7439-96-5	0 - 1.21
Silicon	7440-21-3	0 - 1.21
Titanium dioxide	13463-67-7	0.307 - 1.03
Polybutylate	24936-97-8	0.273 - 0.340
Antimony trioxide	1309-64-4	0.161 - 0.243
Carbon	7440-44-0	0 - 0.201
Boron	7440-42-8	0 - 0.121
Phosphorus	7723-14-0	0 - 0.121
Sulfur	7704-34-9	0 - 0.0806
D&C Green No. 6	128-80-3	0.0394 - 0.0513

Substance	CAS	Model Mass Range (mg)
Silicon dioxide	7631-86-9	0.00422 - 0.00525
Erucamide	112-84-5	0.00110 - 0.00178
4-Dodecylbenzenesulfonic acid	121-65-3	0.000330 - 0.000453

12.0 Summary of Safety and Clinical Performance (SSCP)

The SSCP has been adapted in accordance with the clinical evaluation assessment by the Notified Body on which CE certification has been granted. The SSCP contains a relevant summary of the same information.

The Notified Body has taken notice of and agreed with the benefit-risk rationales for the short- and long-term safety and effectiveness of the SAPIEN 3/SAPIEN 3 Ultra platform.

Conformity with the entire SAPIEN 3/SAPIEN 3 Ultra platform of the Performance Requirements (GSPR) for safety (MDR GSPR1), performance (MDR GSPR1), acceptability of side-effects (MDR GSPR8), usability (MDR GSPR5), device lifetime (MDR GSPR6), and acceptable benefit-risk profile (MDR GSPR8) has been established for the labelled indications.

Refer to <https://meddeviceinfo.edwards.com/> for a SSCP for this medical device.

After the launch of the European Database on Medical Devices/Eudamed, refer to <https://ec.europa.eu/tools/eudamed> for a SSCP for this medical device.

13.0 Basic Unique Device Identification-Device Identifier (UDI-DI)

The Basic UDI-DI is the access key for device-related information entered in the Eudamed. The Basic UDI-DI for the valves, delivery system and sheath can be used to locate the SSCP.

The following tables contain the Basic UDI-DIs:

Table 9: Edwards SAPIEN 3 System

Product	Model				Basic UDI-DI
	20 mm System	23 mm System	26 mm System	29 mm System	
Edwards SAPIEN 3 Transcatheter Heart Valve	9600TFX (20 mm)	9600TFX (23 mm)	9600TFX (26 mm)	9600TFX (29 mm)	0690103D003SAP000VP
Edwards COMMANDER Delivery System	9610TF20	9610TF23	9610TF26	9610TF29	0690103D003COM000TC
Edwards eSheath Introducer Set	9610ES14			9610ES16	0690103D003S3E000NT
Inflation Device	96402			96406	0690103D003IND000TG
Edwards Crimper	9600CR				0690103D003CRI000TH

Table 10: Edwards SAPIEN 3 Ultra System

Product	Model			Basic UDI-DI
	20 mm System	23 mm System	26 mm System	
Edwards SAPIEN 3 Ultra Transcatheter Heart Valve	9750TFX (20 mm)	9750TFX (23 mm)	9750TFX (26 mm)	0690103D00 3SAP000VP
Edwards COMMANDER Delivery System	9610TF20	9610TF23	9610TF26	0690103D00 3COM000TC
Edwards eSheath Introducer Set	9610ES14			0690103D00 3S3E000NT
Inflation Device	96402			0690103D00 3IND000TG
Edwards Crimper	9600CR			0690103D00 3CRI000TH

14.0 Expected Lifetime of the Device

The Edwards Transcatheter Heart Valve has been subjected to rigorous pre-clinical durability testing per the valve testing requirements and in clinical studies and post market studies. The valves were successfully tested to 5 years of simulated wear. In addition, clinical data show durability with follow-up to 5 years. The actual lifetime performance is continuing to be studied and varies from patient to patient.

15.0 Patient Information

A patient implant card is provided with each THV. After implantation, please complete all requested information and provide the implant card to the patient. The serial number is found on the package. This implant card allows patients to inform healthcare providers what type of implant they have when they seek care.

16.0 Recovered THV and Device Disposal

The explanted THV should be placed into a suitable histological fixative such as 10% formalin or 2% glutaraldehyde and returned to the company. Refrigeration is not necessary under these circumstances. Contact Edwards Lifesciences to request an Explant Kit.

Used devices may be handled and disposed of in the same manner as hospital waste and biohazardous materials. There are no special risks related to the disposal of these devices.

17.0 Clinical Studies

Refer to the SSCP for clinical benefits.

18.0 References

1. Bapat V, Attia R, Thomas M. Effect of Valve Design on the Stent Internal Diameter of a Bioprosthetic Valve: A Concept of True Internal Diameter and Its Implications for the Valve-in-Valve Procedure. JACC: Cardiovascular Interventions. Vol. 7, No. 2 2014: 115-127.

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
Figure 1: Edwards SAPIEN 3 Transcatheter Heart Valve

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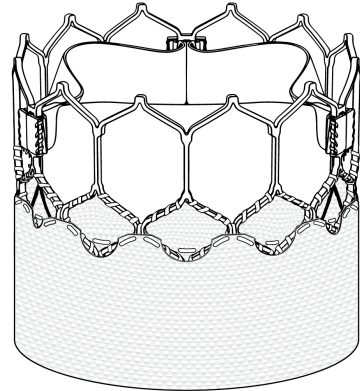Figure 2: Edwards SAPIEN 3 Ultra Transcatheter Heart Valve

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
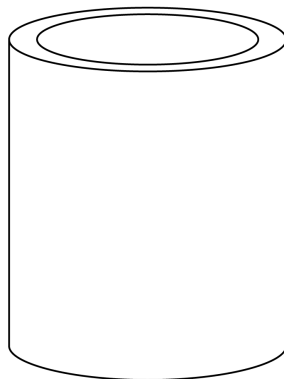
Figure 3: Qualcrimp Crimping Accessory

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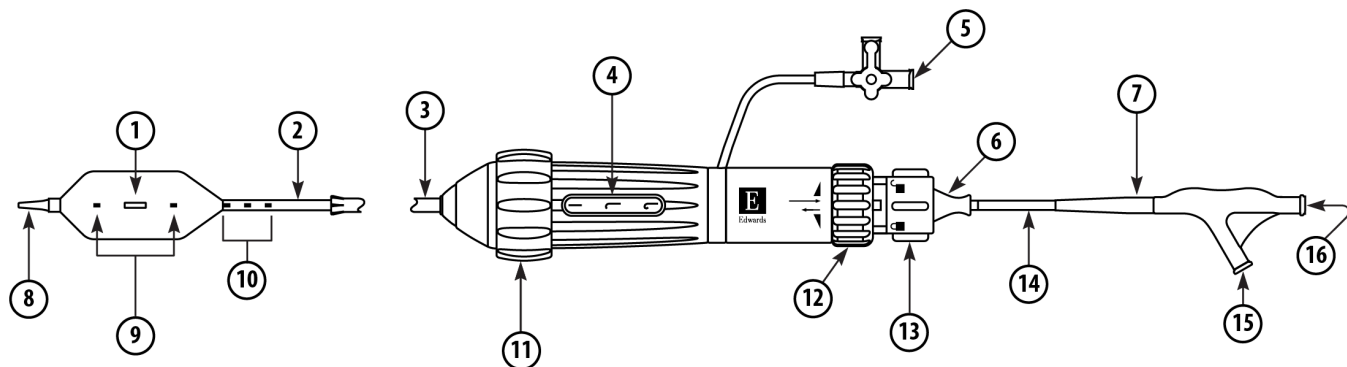	Non-pyrogenic	Ikke-pyrogen	Pyrogeeniton
MD	Medical device	Medisinsk utstyr	Lääkinnällinen laite
	Contains biological material of animal origin	Inneholder biologiske materialer av animalsk opprinnelse	Sisältää eläinperäistä biologista materiaalia
Rx only	Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.	Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.	Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.
	Time & Temperature Sensitive	Tids- og temperaturfølsom	Aika- ja lämpöherkkä
	Contains hazardous substances	Inneholder farlige stoffer	Sisältää vaarallisia aineita
SZ	Size	Størrelse	Koko

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Edwards Lifesciences Services GmbH
Edisonstr. 6
85716 Unterschleissheim
Germany

0344

2022-09
10051155001 A
© Copyright 2022, Edwards Lifesciences LLC
All rights reserved.

Edwards Lifesciences LLC
One Edwards Way
Irvine, CA 92614 USA

Edwards Lifesciences B.V.
Verlengde Poolseweg 16
4818 CL Breda, Netherlands

Telephone +1.949.250.2500
+1.800.424.3278
FAX +1.949.250.2525

Web IFU