Edwards

**Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve System**  
**Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve**  
**Edwards Certitude Delivery System**  
**Edwards SAPIEN 3 og SAPIEN 3 Ultra transkateter hjerteklaffsystem**  
**Edwards SAPIEN 3 og SAPIEN 3 Ultra transkateter hjerteklaff**  
**Edwards Certitude leveringssystem**  
**Katetrilla asennettavat Edwards SAPIEN 3 - ja SAPIEN 3 Ultra -sydänläppäjärjestelmät**  
**Katetrilla asennettavat Edwards SAPIEN 3 - ja SAPIEN 3 Ultra -sydänläpät**  
**Edwards Certitude -asennusjärjestelmä**

| Directory ■ Katalog ■ Hakemisto                             |    |
|-------------------------------------------------------------|----|
| English (en).....                                           | 1  |
| Norsk (no).....                                             | 8  |
| Suomi (fi).....                                             | 15 |
| Figures ■ Figurer ■ Kuvat.....                              | 23 |
| Symbol Legend ■ Symbolforklaring ■ Merkkien selitykset..... | 26 |

## English

### Instructions for Use

Implantation of the transcatheter heart valve should be performed only by physicians who have received Edwards Lifesciences training. The implanting physician should be experienced in standard catheterization techniques. It is at the physician's discretion to choose the appropriate access route to implant the THV based on the patient anatomy and associated risks.

### 1.0 Device Description

#### Edwards SAPIEN 3 and SAPIEN 3 Ultra Transcatheter Heart Valve System

The Edwards SAPIEN 3 and SAPIEN 3 Ultra transcatheter heart valve (THV) system consists of the Edwards SAPIEN 3 and SAPIEN 3 Ultra transcatheter heart valves and delivery systems.

##### • Edwards SAPIEN 3 Transcatheter Heart Valve (Figure 1)

The Edwards SAPIEN 3 transcatheter heart valve (THV) is comprised of a balloon-expandable, radiopaque, cobalt-chromium frame, trileaflet bovine pericardial tissue valve, and polyethylene terephthalate (PET) inner and outer fabric skirts. The leaflets are treated according to the Carpentier-Edwards ThermaFix process.

The THV is intended to be implanted in a native annulus size range associated with the three-dimensional area of the aortic annulus measured at the basal ring during systole as provided in Table 1:

##### • Edwards SAPIEN 3 Ultra Transcatheter Heart Valve (Figure 2)

The Edwards SAPIEN 3 Ultra transcatheter heart valve (THV) is comprised of a balloon-expandable, radiopaque, cobalt-chromium frame, trileaflet bovine pericardial tissue valve, and polyethylene terephthalate (PET) inner and outer fabric skirts. The leaflets are treated according to the Carpentier-Edwards ThermaFix process.

Edwards, Edwards Lifesciences, the stylized E logo, Carpentier-Edwards, Certitude, Edwards SAPIEN, Edwards SAPIEN 3, Edwards SAPIEN 3 Ultra, INSPIRIS, INSPIRIS RESILIA, Qualcrimp, RESILIA, SAPIEN, SAPIEN 3, SAPIEN 3 Ultra, ThermaFix, and VFit are trademarks of Edwards Lifesciences Corporation. All other trademarks are the property of their respective owners.

Table 1

| Native Valve Annulus Size (TEE)* | Native Valve Annulus Size (CT) |                       | THV Size |
|----------------------------------|--------------------------------|-----------------------|----------|
|                                  | Area                           | Area Derived Diameter |          |
| 16 - 19 mm                       | 273 - 345 mm <sup>2</sup>      | 18.6 - 21.0 mm        | 20 mm    |
| 18 - 22 mm                       | 338 - 430 mm <sup>2</sup>      | 20.7 - 23.4 mm        | 23 mm    |
| 21 - 25 mm                       | 430 - 546 mm <sup>2</sup>      | 23.4 - 26.4 mm        | 26 mm    |
| 24 - 28 mm                       | 540 - 683 mm <sup>2</sup>      | 26.2 - 29.5 mm        | 29 mm    |

THV size recommendations are based on native valve annulus size, as measured by transesophageal echocardiography (TEE) or computed tomography (CT). Patient anatomical factors and multiple imaging modalities should be considered during THV size selection.

**Note: Risks associated with undersizing and oversizing should be considered to minimize the risk of paravalvular leak, migration, and/or annular rupture.**

\*Due to limitations in two-dimensional images, 2-D TEE imaging should be supplemented with 3-D area measurements.

Sizing recommendations for implanting the THV in a failing bioprosthesis are provided in Table 2:

Table 2

| Surgical Valve True Inner Diameter (ID) <sup>[1]</sup> | THV Size |
|--------------------------------------------------------|----------|
| 16.5 - 19.0 mm                                         | 20 mm    |
| 18.5 - 22.0 mm                                         | 23 mm    |
| 22.0 - 25.0 mm                                         | 26 mm    |
| 25.0 - 28.5 mm                                         | 29 mm    |

**Note: Surgical valve 'True ID' may be smaller than the labeled valve size. For a failing stentless bioprosthesis, consider sizing recommendations for a native annulus. The dimensions of the failed bioprosthesis should be determined so that the appropriate THV size can be implanted; and is best determined by using computed tomography, magnetic resonance imaging, and/or transesophageal echocardiography.**

Sizing recommendations for implanting the THV in a failing INSPIRIS RESILIA aortic surgical bioprosthesis in sizes 19 - 25 mm, based on bench testing, are provided in Table 3 below:

Table 3

| INSPIRIS RESILIA Aortic Valve (model 11500A)* Labeled Size | THV Size       |
|------------------------------------------------------------|----------------|
| 19 mm                                                      | 20 mm or 23 mm |
| 21 mm                                                      | 23 mm or 26 mm |

| INSPIRIS RESILIA Aortic Valve (model 11500A)* Labeled Size | THV Size |
|------------------------------------------------------------|----------|
| 23 mm                                                      | 26 mm    |
| 25 mm                                                      | 29 mm    |

\*INSPIRIS RESILIA aortic valve model 11500A sizes 19 - 25 mm incorporate VFit technology which consists of expandable bands and fluoroscopically visible size markers designed for potential future valve-in-valve procedures. Clinical data are not currently available on the INSPIRIS RESILIA aortic valve Model 11500A valve-in-valve procedure or expansion feature. The impact of tissue ingrowth on the expansion feature of the INSPIRIS RESILIA aortic valve has not been assessed.

**WARNING: Do not perform stand-alone balloon aortic valvuloplasty procedures in the INSPIRIS RESILIA aortic valve for the sizes 19 - 25 mm. This may expand the valve causing aortic incompetence, coronary embolism or annular rupture.**

**Note: INSPIRIS RESILIA aortic valve model 11500A sizes 27 - 29 mm do not incorporate VFit technology and therefore follow the surgical valve True ID sizing provided in Table 2.**

**Note: Exact volume required to deploy the THV may vary depending on the bioprosthesis inner diameter. Factors such as calcification and pannus tissue growth may not be accurately visualized in imaging and may reduce the effective inner diameter of the failing bioprosthesis to a size smaller than the 'True ID'. These factors should be considered and assessed in order to determine the most appropriate THV size to achieve nominal THV deployment and sufficient anchoring. Do not exceed the rated burst pressure. See inflation parameters in Table 4.**

#### • Edwards Certitude Delivery System (Figure 3)

The Edwards Certitude delivery system facilitates the placement of the bioprosthesis. The delivery system consists of a flex catheter to aid in tracking and valve positioning. The delivery system includes a tapered tip to facilitate crossing of the valve. The handle contains a flex wheel to control flexing of the balloon catheter. A stylet is included within the guidewire lumen of the delivery system. A radiopaque center marker in the balloon is provided to help with valve positioning. The extension tubing is used during THV deployment.

The inflation parameters for valve deployment are:

**Table 4**

| Model    | Nominal Balloon Diameter | Nominal Inflation Volume | Rated Burst Pressure (RBP) |
|----------|--------------------------|--------------------------|----------------------------|
| 9620TA20 | 20 mm                    | 12 ml                    | 7 atm                      |
| 9620TA23 | 23 mm                    | 17 ml                    | 7 atm                      |
| 9620TA26 | 26 mm                    | 23 ml                    | 7 atm                      |
| 9620TA29 | 29 mm                    | 30 ml                    | 7 atm                      |

#### • Edwards Certitude Introducer Sheath Set (Figure 4)

The Edwards Certitude introducer sheath set facilitates the introduction and removal of devices utilized with the SAPIEN 3 and SAPIEN 3 Ultra transcatheter heart valve. The sheath has a radiopaque marker for visualization of the sheath tip and non-radiopaque depth markings on the distal end of the body of the sheath. The proximal end of the sheath includes a flush tube and three hemostasis valves. An introducer is supplied with the sheath. The entire introducer is radiopaque.

**Table 5**

#### Introducer Sheath Set Information

| Model                                              | 9620IS18           | 9620IS21     |
|----------------------------------------------------|--------------------|--------------|
| Sheath Inside Diameter                             | 18F (6.1 mm)       | 21F (6.9 mm) |
| Sheath Effective Length                            | 21 cm              | 21 cm        |
| Introducer Size                                    | OD: 6.3 mm         | OD: 7.0 mm   |
| Introducer Effective Length                        | 33 cm              |              |
| Diameter of the largest guidewire that can be used | 0.035 in (0.89 mm) |              |

#### • Qualcrimp Crimping Accessory

The Qualcrimp crimping accessory is used during THV crimping (Figure 5).

#### • Loader

The loader is used to aid insertion of the delivery system into the sheath (Figure 6).

#### • Edwards Crimper and Crimp Stopper (Figure 7)

The Edwards crimper reduces the diameter of the valve to mount it onto its delivery system. The crimper is comprised of a housing and a compression mechanism that is closed with a handle located on the housing. A 2-piece crimp stopper is used to crimp the valve to its intended diameter.

#### • Inflation Device

An inflation device with locking mechanism is used during valve deployment.

**Note: For proper volume sizing, the delivery system must be used with the inflation device provided by Edwards Lifesciences.**

## 2.0 Intended Use

The bioprosthesis is intended for use in patients requiring heart valve replacement. The delivery system and accessories are intended to facilitate the placement of the bioprosthesis via the transapical and transaortic access approaches.

## 3.0 Indications

1. The Edwards SAPIEN 3 Ultra transcatheter heart valve system is indicated for use in patients with severe, symptomatic, calcific aortic valve stenosis who are judged by a Heart Team to be at intermediate or greater risk for open heart surgery (i.e., predicted risk of surgical mortality  $\geq 3\%$  at 30 days, based on the Society of Thoracic Surgeons (STS) risk score and other clinical comorbidities unmeasured by the STS risk calculator).
2. The Edwards SAPIEN 3 Ultra transcatheter heart valve system is indicated for patients with symptomatic heart disease due to failure (stenosed, insufficient, or combined) of a surgical bioprosthetic aortic valve or a surgical bioprosthetic mitral valve who are judged by a heart team, including a cardiac surgeon, to be at high or greater risk for open surgical therapy (i.e., predicted risk of surgical mortality  $\geq 8\%$  at 30 days, based on the Society of Thoracic Surgeons (STS) risk score and other clinical co-morbidities unmeasured by the STS risk calculator).

## 4.0 Contraindications

Use of the Edwards SAPIEN 3 and SAPIEN 3 Ultra THV system is contraindicated in patients who:

- Cannot tolerate anticoagulation/antiplatelet regimen or who have active bacterial endocarditis or other active infections.

## 5.0 Warnings

- The devices are designed, intended, and distributed STERILE for single use only. **Do not resterilize or reuse the devices.** There are no data to support the sterility, nonpyrogenicity, and functionality of the devices after reprocessing.
- Correct sizing of the THV is essential to minimize the risk of paravalvular leak, migration, and/or annular rupture.
- The physician must verify correct orientation of the THV prior to its implantation.
- Accelerated deterioration of the THV may occur in patients with an altered calcium metabolism.
- When using venous pacing, observation of the pacing lead throughout the procedure is essential to avoid the potential risk of the pacing lead causing a cardiovascular perforation.
- The THV must remain hydrated at all times and cannot be exposed to solutions, antibiotics, chemicals, etc. other than its shipping storage solution and sterile physiologic saline solution to prevent leaflet damage that may impact valve functionality. THV leaflets mishandled or damaged during any part of the procedure will require replacement of the THV.
- Patients with hypersensitivities to cobalt, nickel, chromium, molybdenum, titanium, manganese, silicon, bovine tissue, and/or polymeric materials may have an allergic reaction to these materials.
- Do not use the THV if the tamper evident seal is broken, as sterility may be compromised.

- Do not use the THV if the temperature indicator has been activated, as valve function may be compromised.
- Do not use the THV if the expiration date has elapsed, as either sterility or valve function may be compromised.
- Do not use the THV if the storage solution does not completely cover the valve or the valve is damaged.
- Do not mishandle the delivery system or use the delivery system and accessory devices if the packaging sterile barriers and any components have been opened or damaged, cannot be flushed, or the expiration date has elapsed.
- Valve recipients should be maintained on anticoagulant/antiplatelet therapy, except when contraindicated, to minimize the risk of valve thrombosis or thromboembolic events, as determined by their physicians. This device has not been tested for use without anticoagulation.
- The procedure should be conducted under fluoroscopic guidance. Some fluoroscopically guided procedures are associated with a risk of radiation injury to the skin. These injuries may be painful, disfiguring, and long-lasting.
- Caution should be exercised in implanting a valve in patients with clinically significant coronary artery disease.
- Patients with pre-existing bioprostheses should be carefully assessed prior to implantation of the valve to ensure proper valve positioning and deployment.
- Balloon valvuloplasty should be avoided in the treatment of failing bioprostheses as this may result in embolization of bioprosthesis material and mechanical disruption of the valve leaflets.

## 6.0 Precautions

- Long-term durability has not been established for the THV. Regular medical follow-up is advised to evaluate valve performance.
- Glutaraldehyde may cause irritation of the skin, eyes, nose and throat. Avoid prolonged or repeated exposure to, or breathing of, the solution. Use only with adequate ventilation. If skin contact occurs, immediately flush the affected area with water; in the event of contact with eyes, seek immediate medical attention. For more information about glutaraldehyde exposure, refer to the Material Safety Data Sheet available from Edwards Lifesciences.
- The safety and effectiveness of the THV implantation has not been established in patients who have:
  - Non-calcified aortic annulus
  - Congenital unicuspid aortic valve
  - Pre-existing prosthetic ring in any position
  - Severe ventricular dysfunction with ejection fraction < 20%
  - Hypertrophic cardiomyopathy with or without obstruction
  - Aortic stenosis characterized by a combination of AV low flow, low gradient
  - A concomitant paravalvular leak where the failing bioprosthesis is not securely fixed in the native annulus or is not structurally intact (e.g., wireframe fracture)
  - A partially detached leaflet of the failing bioprosthesis that in the aortic position may obstruct a coronary ostium
  - Bulky calcified aortic valve leaflets in close proximity to coronary ostia.
- Appropriate antibiotic prophylaxis is recommended post-procedure in patients at risk for prosthetic valve infection and endocarditis.
- Special care must be exercised in mitral valve replacement if chordal preservation techniques were used in the primary implantation to avoid entrapment of the subvalvular apparatus.
- Based on the treating physician's consideration of risks and benefits, the valve may be implanted in relatively young patients, although the longer-term durability is still the subject of ongoing clinical research.
- Do not overinflate the deployment balloon, as this may prevent proper valve leaflet coaptation and thus impact valve functionality.
- THV recipients should be maintained on anticoagulant/antiplatelet therapy to minimize risk of valve thrombosis or thromboembolic events, as determined by their physicians.
- Residual mean gradient may be higher in a "THV-in-failing bioprosthesis" configuration than that observed following implantation of the valve inside a native aortic annulus using the same size device. Patients with elevated mean gradient post procedure should be carefully followed. It is important that the manufacturer, model and size of the preexisting bioprosthetic valve be determined, so that the appropriate valve can be implanted and a prosthesis-patient mismatch be avoided. Additionally, pre-procedure imaging

modalities must be employed to make as accurate a determination of the inner diameter as possible.

## 7.0 Potential Adverse Events

Potential risks associated with the overall procedure including access, cardiac catheterization, and general anesthesia:

- Death
- Stroke/transient ischemic attack, clusters or neurological deficit
- Paralysis
- Permanent disability
- Respiratory insufficiency or respiratory failure
- Cardiovascular injury including perforation or dissection of vessels, ventricle, atrium, septum, myocardium or valvular structures that may require intervention
- Pericardial effusion or cardiac tamponade
- Thoracic bleeding
- Embolization including air, calcific valve material or thrombus
- Infection including septicemia and endocarditis
- Heart failure
- Myocardial ischemia or infarction
- Renal insufficiency or renal failure
- Conduction system defect which may require a permanent pacemaker
- Arrhythmias including ventricular fibrillation (VF) and ventricular tachycardia (VT)
- Retroperitoneal bleed
- Arteriovenous (AV) fistula or pseudoaneurysm
- Reoperation
- Ischemia or nerve injury or brachial plexus injury or compartment syndrome
- Restenosis
- Pulmonary edema
- Pleural effusion
- Bleeding, bleeding requiring transfusion or intervention
- Anemia
- Vessel thrombosis/occlusion
- Abnormal lab values (including electrolyte imbalance)
- Hypertension or hypotension
- Allergic reaction to anesthesia, contrast media, device materials or bovine pericardial tissue
- Hematoma
- Syncope
- Pain or changes (e.g., wound infection, hematoma, and other wound care complications) at the access site
- Exercise intolerance or weakness
- Inflammation
- Angina
- Vasovagal response
- Heart murmur
- Fever

Additional potential risks associated with the TAVR procedure, the bioprosthesis, and the use of its associated devices and accessories include:

- Cardiac arrest
- Cardiogenic shock
- Emergency cardiac surgery
- Cardiac failure or low cardiac output
- Coronary flow obstruction/transvalvular flow disturbance
- Device thrombosis requiring intervention
- Valve thrombosis
- Device embolization
- Device migration or malposition requiring intervention
- Left ventricular outflow tract obstruction
- Valve deployment in unintended location
- Valve stenosis
- Structural valve deterioration (wear, fracture, calcification, leaflet tear/tearing from the stent posts, leaflet retraction, suture line disruption of components of a prosthetic valve, thickening, stenosis)
- Device degeneration
- Paravalvular or transvalvular leak
- Valve regurgitation

- Hemolysis
- Device explants
- Mediastinitis
- Mediastinal bleeding
- Nonstructural dysfunction
- Mechanical failure of delivery system, and/or accessories, including balloon rupture and tip separation
- Non-emergent reoperation
- Allergic/immunologic reaction to the implant
- Injury to the mitral valve

For a patient/user/third party in the European Economic area; if, during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and your national competent authority, which can be found at [https://ec.europa.eu/growth/sectors/medical-devices/contacts\\_en](https://ec.europa.eu/growth/sectors/medical-devices/contacts_en).

## 8.0 Directions for Use

### 8.1 System Compatibility

Table 6

|                                                                                                           | 20 mm<br>System | 23 mm<br>System | 26 mm<br>System | 29 mm<br>System |
|-----------------------------------------------------------------------------------------------------------|-----------------|-----------------|-----------------|-----------------|
| Product Name                                                                                              | Model           |                 |                 |                 |
| Edwards SAPIEN 3 transcatheter heart valve                                                                | 9600TFX (20 mm) | 9600TFX (23 mm) | 9600TFX (26 mm) | 9600TFX (29 mm) |
| Edwards Certitude delivery system                                                                         | 9620TA20        | 9620TA23        | 9620TA26        | 9620TA29        |
| Edwards Certitude introducer sheath set                                                                   | 9620IS18 (18F)  |                 |                 | 9620IS21 (21F)  |
| inflation device                                                                                          | 96402           |                 |                 | 96406           |
| Edwards crimper                                                                                           | 9600CR          |                 |                 |                 |
| Qualcrimp crimping accessory, crimp stopper, loader and extension tubing provided by Edwards Lifesciences |                 |                 |                 |                 |

Table 7

|                                                                                                              | 20 mm<br>System    | 23 mm<br>System    | 26 mm<br>System    |
|--------------------------------------------------------------------------------------------------------------|--------------------|--------------------|--------------------|
| Product Name                                                                                                 | Model              |                    |                    |
| Edwards SAPIEN 3<br>Ultra transcatheter<br>heart valve                                                       | 9750TFX<br>(20 mm) | 9750TFX<br>(23 mm) | 9750TFX<br>(26 mm) |
| Edwards Certitude<br>delivery system                                                                         | 9620TA20           | 9630TA23           | 9630TA26           |
| Edwards Certitude<br>introducer sheath set                                                                   | 9620IS18           | 9620IS21           |                    |
| inflation device                                                                                             | 96402              |                    |                    |
| Edwards crimper                                                                                              | 9600CR             |                    |                    |
| Qualcrimp crimping accessory, crimp stopper, loader and extension<br>tubing provided by Edwards Lifesciences |                    |                    |                    |

#### Additional Equipment

- Balloon catheter per the discretion of the physician
- 20 cc syringe or larger
- 50 cc syringe or larger
- High-pressure 3-way stopcock
- Standard cardiac catheterization lab equipment and supplies, and access to standard heart valve operating room equipment and supplies
- Fluoroscopy (fixed, mobile or semi-mobile fluoroscopy systems appropriate for use in percutaneous coronary interventions)
- Transesophageal or transthoracic echocardiography capabilities

- 18 gauge Seldinger needle (for transaortic)
- 145 cm x 0.035 in (0.89 mm) soft guidewire
- 180 cm or 260 cm x 0.035 in (0.89 mm) & Exchange length 0.035 in (0.89 mm) extra-stiff guidewires
- Temporary pacemaker (PM) and pacing lead
- Sterile rinsing basins, physiological saline, heparinized saline, 15% diluted radiopaque contrast medium
- Sterile table for THV and accessories preparation

### 8.2 THV Handling and Preparation

Follow sterile technique during device preparation and implantation.

#### 8.2.1 THV Rinsing Procedure

Before opening the valve jar, carefully examine for evidence of damage (e.g., a cracked jar or lid, leakage, or broken or missing seals).

**CAUTION: If the container is found to be damaged, leaking, without adequate sterilant, or missing intact seals, the THV must not be used for implantation, as sterility may be compromised.**

- Set up two (2) sterile bowls with at least 500 ml of sterile physiological saline to thoroughly rinse the THV.
- Carefully remove the valve/holder assembly from the jar without touching the tissue. Verify the valve serial identification number with the number on the jar lid and record in the patient information documents. Inspect the valve for any signs of damage to the frame or tissue.
- Rinse the THV as follows:
  - Place the THV in the first bowl of sterile, physiological saline. Be sure the saline solution completely covers the THV and holder.
  - With the valve and holder submerged, slowly agitate (to gently swirl the valve and holder) back and forth for a minimum of 1 minute.
  - Transfer the THV and holder to the second rinsing bowl of sterile physiological saline and gently agitate for at least one more minute. Ensure the rinse solution in the first bowl is not used.
  - The valve should be left in the final rinse solution until needed to prevent the tissue from drying.

**CAUTION: Do not allow the valve to come into contact with the bottom or sides of the rinse bowl during agitation or swirling in the rinse solution. Direct contact between the identification tag and valve is also to be avoided during the rinse procedure. No other objects should be placed in the rinse bowls. The valve should be kept hydrated to prevent the tissue from drying.**

#### 8.2.2 Prepare the System

- Visually inspect all components for damage. Ensure the system is fully unflexed.
- Prime and flush the introducer and sheath with heparinized saline. Hydrate the length of the introducer and sheath.
- Advance the introducer fully into the sheath housing.
- Unscrew the loader cap from the loader and flush the loader cap with heparinized saline.
- Place the loader cap onto the delivery system with the inside of the cap oriented towards the tapered tip.
- Flush the extension tubing and connect to the delivery system.
- Partially fill a 50 ml or larger syringe with diluted contrast medium, and connect to the extension tubing.
- Fill the inflation device with 20 ml of diluted contrast medium, lock the inflation device, and connect to the extension tubing. Close 3-way stopcock to inflation device.
- De-air the delivery system using the luer lock syringe. Leave zero-pressure in the system. Close the 3-way stopcock to the luer lock syringe.
- Remove 3 ml fluid from the delivery system by turning the knob of the locked inflation device. Keep the inflation device locked for THV crimping steps.

#### 8.2.3 Mount and Crimp the THV onto the Delivery System

- Completely submerge the Qualcrimp crimping accessory in a bowl of 100 ml physiological saline solution. Gently compress until fully saturated. Swirl for a minimum of 1 minute. Repeat this process in a second bowl.
- Rotate the crimper until the aperture is fully opened. Attach the 2-piece crimp stopper to the crimper.


3. Remove the THV from the holder and remove the ID tag.
4. If necessary, partially crimp the THV in the crimper until it fits snugly inside the Qualcrimp crimping accessory.

**Note: Partial crimping is not necessary for the 20 mm valve.**

5. Place the Qualcrimp crimping accessory over the THV.
6. The orientation of the THV on the delivery system is described below:


#### Antegrade Approach:

Inflow (outer skirt end) of the THV towards the **proximal end** of the delivery system.
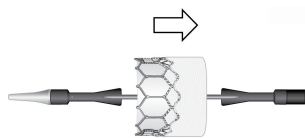
#### Retrograde Approach:

Inflow (outer skirt end) of the THV towards the **distal end** of the delivery system.

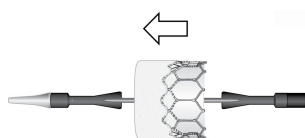7. Place the THV and Qualcrimp crimping accessory in the crimper. Insert the delivery system coaxially into the THV.
  8. Crimp the THV between the two internal shoulders of the delivery system until it reaches the Qualcrimp stop.
  9. Remove the Qualcrimp crimping accessory from the THV/balloon assembly and Qualcrimp stop from the crimp stopper, leaving the final stop in place.
- Note: Ensure that the THV remains centered and coaxial within the two internal shoulders.**
10. Place the THV/balloon assembly back in the crimper aperture, fully crimp the THV until it reaches the final stop and hold for 5 seconds.
  11. Repeat the full crimp of the THV two times for a total of 3 crimps.
  12. Flush the loader with heparinized saline. Immediately advance the loader over the THV until the tapered tip of the delivery system is exposed and the THV is within the distal end of the loader tube.

**CAUTION: The THV should not remain fully crimped and/or in the loader over 15 minutes, as leaflet damage may result and impact valve functionality.**

13. Attach the loader cap to the loader and flush through the flush port on the loader. Remove the stylet and flush the guidewire lumen of the delivery system.

**CAUTION: Keep the THV hydrated until ready for implantation to prevent damage to the leaflets which may impact valve functionality.**

**WARNING: The physician must verify correct orientation of the THV prior to its implantation to prevent the risk of severe patient harm.**

14. With 3-way stopcock still closed to the luer lock syringe, unlock the inflation device. Allow the delivery system to reach zero-pressure.
15. Close the 3-way stopcock to the delivery system. Use the luer lock syringe to de-air the inflation device if necessary.
16. Adjust the inflation device to the inflation volume required to deploy the THV, per table 4.  
Re-lock the inflation device. Close the 3-way stopcock to the luer lock syringe and remove syringe.

**CAUTION: Maintain the inflation device in a locked position until THV deployment to prevent premature balloon inflation and subsequent improper THV deployment.**

### 8.3 Native Valve Predilation and THV Delivery

Native valve predilation and THV delivery should be performed under general anesthesia with hemodynamic monitoring in a catheterization

lab/hybrid operating room with fluoroscopic and echocardiographic imaging capabilities.

The following tables show the minimum required distances from the valvular plane to the distal tip of the sheath to allow the Edwards Certitude delivery system balloon to inflate properly during THV deployment. **These distances do not include sheath insertion depth**, which should be considered during the transaortic approach when selecting the access site on the ascending aorta.

**Table 8: Edwards SAPIEN 3 System**

| Delivery System | THV   | Minimum Required Distance From Sheath Tip to Valvular Plane |
|-----------------|-------|-------------------------------------------------------------|
| Model 9620TA20  | 20 mm | 3.5 cm                                                      |
| Model 9620TA23  | 23 mm | 3.5 cm                                                      |
| Model 9620TA26  | 26 mm | 3.5 cm                                                      |
| Model 9620TA29  | 29 mm | 4.0 cm                                                      |

**Table 9: Edwards SAPIEN 3 Ultra System**

| Delivery System | THV   | Minimum Required Distance From Sheath Tip to Valvular Plane |
|-----------------|-------|-------------------------------------------------------------|
| Model 9620TA20  | 20 mm | 3.5 cm                                                      |
| Model 9630TA23  | 23 mm | 3.5 cm                                                      |
| Model 9630TA26  | 26 mm | 3.5 cm                                                      |

Administer heparin to maintain the ACT at  $\geq 250$  sec.

**CAUTION: Contrast media use should be monitored to reduce the risk of renal injury.**

#### 8.3.1 Baseline Parameters

1. Advance a 5F (1.67 mm) or 6F (2.0 mm) pigtail catheter into the descending aorta and perform a supra-aortic angiogram with the projection of the native aortic valve perpendicular to the view.
2. Evaluate the distances of the right and left coronary ostia from the aortic annulus in relation to the THV frame height.
3. Introduce a pacemaker (PM) lead and position appropriately.
4. Set the stimulation parameters to obtain 1:1 capture, and test pacing.

#### 8.3.2 Access

**CAUTION: Care should be taken to avoid damage to soft tissue, chordae, aorta, native leaflet or ventricular wall during insertion, positioning and removal of devices.**

##### 8.3.2.1 Transapical Access

1. Access the apex through an anterior mini thoracotomy at the 5th or 6th intercostal space. Incise the pericardium to expose the apex of the left ventricle (LV).
2. Attach epicardial pacing leads to left ventricle or insert transvenous pacing leads and secure proximal ends of leads into pacemaker. Set the stimulation parameters, test rapid pacing.
3. Place a reinforced double purse string on the LV apex to access the left ventricle.
4. Gain aortic valve access through standard transapical techniques.
5. Using the sheath depth markers, advance the introducer and sheath over the guidewire to the desired depth while following its progression on fluoroscopy.
6. Withdraw the introducer slowly, keeping the sheath in place. Maintain guidewire position across the valve.

##### 8.3.2.2 Transaortic Access

1. Access the ascending aorta using standard surgical technique (e.g., a partial J-sternotomy or right parasternal mini thoracotomy).
2. Place two reinforced purse string sutures at the intended access site in the ascending aorta.

**Note: The selected access site should be soft by digital palpation.**

3. Introduce a pacemaker lead until its distal end is positioned in the right ventricle. Set the stimulation parameters and test pacing.
4. Gain aortic valve access through standard transaortic techniques.
5. Insert the Edwards Certitude introducer sheath set, or desired introducer sheath for BAV, into the aorta to approximately 2 cm. Withdraw the introducer slowly, keeping the sheath in place. Maintain guidewire position across the aortic valve.

### 8.3.3 Valvuloplasty

Pre-dilate the native aortic valve, per the discretion of the physician, according to the instructions for use for the selected balloon aortic valvuloplasty catheter.

**CAUTION: Valve implantation should not be carried out if the balloon cannot be fully inflated during valvuloplasty.**

### 8.3.4 THV Delivery

**CAUTION: Care should be taken to avoid damage to soft tissue, chordae, aorta, native leaflet or ventricular wall during insertion, positioning and removal of devices.**

1. Confirm that the THV is oriented properly and the volume in the inflation device matches the indicated volume.
2. Advance the THV/balloon assembly with the loader over the guidewire.
3. Engage loader into the sheath housing while maintaining a firm grip.
4. Advance the valve out of the loader into the large section of the sheath. Tap on the sheath housing to release air bubbles to the proximal end of the loader. Depress button valve on loader to de-air.
5. Advance the THV/balloon assembly through the sheath and position within the target valve. If needed, rotate the flex wheel on the handle to articulate the THV/balloon assembly into position.

**CAUTION: To prevent possible leaflet damage that may impact valve functionality, the THV should not remain in the sheath for over 5 minutes.**

6. Ensure that the THV is correctly positioned between the two internal shoulders of the delivery system.
7. Begin THV deployment:
  - a) Unlock the inflation device provided by Edwards Lifesciences.
  - b) Begin rapid pacing; once systolic blood pressure has decreased to 50 mmHg or below, balloon inflation can commence.
  - c) Deploy the valve by inflating the balloon with the entire volume in the Inflation device provided by Edwards Lifesciences, hold for 3 seconds and confirm that the barrel of the inflation device is empty to ensure complete inflation of the balloon.
  - d) Deflate the balloon. When the balloon has been completely deflated, turn off the pacemaker.

### 8.3.5 System Removal

1. If articulation was used, completely unflex the delivery system. Retract the delivery system and guidewire into the sheath. Remove the loader and delivery system from the sheath.

**CAUTION: Properly deflate the balloon and unflex the delivery system prior to removal.**

2. Remove all devices when the ACT level is appropriate.
3. Remove the sheath from the access site, close the access site and confirm hemostasis.

### 8.4 Verification of THV Position and Measurements

Measure and record hemodynamic parameters.

1. Perform an angiogram to evaluate device performance and coronary patency, where applicable.
2. Measure and record the transvalvular pressure gradients.
3. Remove all devices when the ACT level is appropriate (e.g., reaches < 150 sec). Refer to the introducer sheath instructions for use for device removal.
4. Close the access site.

## 9.0 How Supplied

STERILE: The valve is supplied sterilized with glutaraldehyde solution.

The delivery system and accessories are supplied sterilized by ethylene oxide gas.

The THV is supplied nonpyrogenic packaged in buffered glutaraldehyde, in a plastic jar to which a tamper evident seal has been applied. Each jar is shipped in a shelf box containing a temperature indicator to detect exposure of the THV to extreme temperature. The shelf box is enclosed in Styrofoam prior to shipping.

### 9.1 Storage

The THV must be stored at 10 °C to 25 °C (50 °F to 77 °F). Each jar is shipped in an enclosure containing a temperature indicator to detect exposure of the THV to extreme temperature.

The delivery system and accessories should be stored in a cool, dry place.

## 10.0 MR Safety

MR Conditional

Non-clinical testing has demonstrated that the Edwards SAPIEN 3 Ultra transcatheter heart valve and Edwards SAPIEN 3 transcatheter heart valve are MR Conditional. A patient with this device can be scanned safely, immediately after placement of this device under the following conditions:

- Static magnetic field of 1.5 tesla (T) or 3.0 tesla (T)
- Maximum spatial gradient field of 2500 Gauss/cm (25 T/m) or less
- Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2.0 W/kg (Normal Operating Mode)

Under the scan conditions defined above, the transcatheter heart valve is expected to produce a maximum temperature rise of 3.0 °C after 15 minutes of continuous scanning.

In non-clinical testing, the image artifact caused by the device extends as far as 14.5 mm from the implant for spin echo images and 30 mm for gradient echo images when scanned in a 3.0 T MRI system. The artifact obscures the device lumen in gradient echo images.

The implant has not been evaluated in MR systems other than 1.5 T or 3.0 T.

For valve-in-valve implantation or in the presence of other implants, please refer to the MRI safety information for the surgical valve or other devices prior to MR imaging.

## 11.0 Qualitative and Quantitative Information related to the THV

This device contains the following substance(s) defined as CMR 1B in a concentration above 0.1% weight by weight:

Cobalt; CAS No. 7440-48-4; EC No. 231-158-0

Current scientific evidence supports that medical devices manufactured from cobalt alloys or stainless steel alloys containing cobalt do not cause an increased risk of cancer or adverse reproductive effects.

For THV, the following tables show the qualitative and quantitative information on the materials and substances:

Table 10: SAPIEN 3 Transcatheter Heart Valve

| Substance                                       | CAS          | Model Mass Range (mg) |
|-------------------------------------------------|--------------|-----------------------|
| Cobalt                                          | 7440-48-4    | 131 - 427             |
| Nickel                                          | 7440-02-0    | 148 - 405             |
| Chromium                                        | 7440-47-3    | 85.2 - 230            |
| Polyethylene terephthalate                      | 25038-59-9   | 102 - 170             |
| Collagens, bovine, polymers with glutaraldehyde | 2370819-60-4 | 58.3 - 141            |
| Molybdenum                                      | 7439-98-7    | 40.3 - 115            |
| Polytetrafluoroethylene                         | 9002-84-0    | 17.5 - 25.5           |
| Polyethylene                                    | 9002-88-4    | 14.2 - 19.7           |
| Iron                                            | 7439-89-6    | 0 - 10.9              |
| Titanium                                        | 7440-32-6    | 0 - 10.9              |
| Manganese                                       | 7439-96-5    | 0 - 1.64              |

| Substance                     | CAS        | Model Mass Range (mg) |
|-------------------------------|------------|-----------------------|
| Silicon                       | 7440-21-3  | 0 - 1.64              |
| Titanium dioxide              | 13463-67-7 | 0.219 - 0.752         |
| Polybutylate                  | 24936-97-8 | 0.273 - 0.383         |
| Carbon                        | 7440-44-0  | 0 - 0.274             |
| Antimony trioxide             | 1309-64-4  | 0.112 - 0.190         |
| Boron                         | 7440-42-8  | 0 - 0.164             |
| Phosphorus                    | 7723-14-0  | 0 - 0.164             |
| Sulfur                        | 7704-34-9  | 0 - 0.109             |
| D&C Green No. 6               | 128-80-3   | 0.0394 - 0.0578       |
| Silicon dioxide               | 7631-86-9  | 0.00422 - 0.00592     |
| Erucamide                     | 112-84-5   | 0.000683 - 0.00128    |
| 4-Dodecylbenzenesulfonic acid | 121-65-3   | 0.000286 - 0.000430   |

**Table 11: SAPIEN 3 Ultra Transcatheter Heart Valve**

| Substance                                       | CAS          | Model Mass Range (mg) |
|-------------------------------------------------|--------------|-----------------------|
| Cobalt                                          | 7440-48-4    | 131 - 314             |
| Nickel                                          | 7440-02-0    | 148 - 298             |
| Polyethylene terephthalate                      | 25038-59-9   | 142 - 212             |
| Chromium                                        | 7440-47-3    | 85.2 - 169            |
| Collagens, bovine, polymers with glutaraldehyde | 2370819-60-4 | 58.3 - 97.5           |
| Molybdenum                                      | 7439-98-7    | 40.3 - 84.6           |
| Polyethylene                                    | 9002-88-4    | 19.4 - 22.0           |
| Polytetrafluoroethylene                         | 9002-84-0    | 12.3 - 15.1           |
| Iron                                            | 7439-89-6    | 0 - 8.06              |
| Titanium                                        | 7440-32-6    | 0 - 8.06              |
| Manganese                                       | 7439-96-5    | 0 - 1.21              |
| Silicon                                         | 7440-21-3    | 0 - 1.21              |
| Titanium dioxide                                | 13463-67-7   | 0.307 - 1.03          |
| Polybutylate                                    | 24936-97-8   | 0.273 - 0.340         |
| Antimony trioxide                               | 1309-64-4    | 0.161 - 0.243         |
| Carbon                                          | 7440-44-0    | 0 - 0.201             |
| Boron                                           | 7440-42-8    | 0 - 0.121             |
| Phosphorus                                      | 7723-14-0    | 0 - 0.121             |
| Sulfur                                          | 7704-34-9    | 0 - 0.0806            |
| D&C Green No. 6                                 | 128-80-3     | 0.0394 - 0.0513       |
| Silicon dioxide                                 | 7631-86-9    | 0.00422 - 0.00525     |
| Erucamide                                       | 112-84-5     | 0.00110 - 0.00178     |
| 4-Dodecylbenzenesulfonic acid                   | 121-65-3     | 0.000330 - 0.000453   |

## 12.0 Summary of Safety and Clinical Performance (SSCP)

The SSCP has been adapted in accordance with the clinical evaluation assessment by the Notified Body on which CE certification has been granted. The SSCP contains a relevant summary of the same information.

The Notified Body has taken notice of and agreed with the benefit-risk rationales for the short- and long-term safety and effectiveness of the SAPIEN 3/SAPIEN 3 Ultra platform.

Conformity with the entire SAPIEN 3/SAPIEN 3 Ultra platform of the performance requirements (GSPR) for safety (MDR GSPR1), performance (MDR GSPR1), acceptability of side-effects (MDR GSPR8), usability (MDR

GSPR5), device lifetime (MDR GSPR6), and acceptable benefit-risk profile (MDR GSPR8) has been established for the labelled indications.

Refer to <https://meddeviceinfo.edwards.com/> for a SSCP for this medical device.

After the launch of the European Database on Medical Devices/Eudamed, refer to <https://ec.europa.eu/tools/eudamed> for a SSCP for this medical device.

## 13.0 Basic Unique Device Identification-Device Identifier (UDI-DI)

The Basic UDI-DI is the access key for device-related information entered in the Eudamed. The Basic UDI-DI for the valves, delivery system and sheath can be used to locate the SSCP.

The following table contains the Basic UDI-DI:

**Table 12: Edwards SAPIEN 3 System**

| Product                                    | Model           |                 |                 |                 | Basic UDI-DI        |
|--------------------------------------------|-----------------|-----------------|-----------------|-----------------|---------------------|
|                                            | 20 mm System    | 23 mm System    | 26 mm System    | 29 mm System    |                     |
| Edwards SAPIEN 3 transcatheter heart valve | 9600TFX (20 mm) | 9600TFX (23 mm) | 9600TFX (26 mm) | 9600TFX (29 mm) | 0690103D003SAP000VP |
| Edwards Certitude delivery system          | 9620TA20        | 9620TA23        | 9620TA26        | 9620TA29        | 0690103D003CER000QZ |
| Edwards Certitude introducer sheath set    | 9620IS18        |                 |                 | 9620IS21        | 0690103D003CIS000SL |
| inflation device                           | 96402           |                 |                 | 96406           | 0690103D003IND000TG |
| Edwards crimper                            | 9600CR          |                 |                 |                 | 0690103D003CRI000TH |

**Table 13: Edwards SAPIEN 3 Ultra System**

| Product                                          | Model           |                 |                 | Basic UDI-DI        |
|--------------------------------------------------|-----------------|-----------------|-----------------|---------------------|
|                                                  | 20 mm System    | 23 mm System    | 26 mm System    |                     |
| Edwards SAPIEN 3 Ultra transcatheter heart valve | 9750TFX (20 mm) | 9750TFX (23 mm) | 9750TFX (26 mm) | 0690103D003SAP000VP |
| Edwards Certitude delivery system                | 9620TA20        | 9630TA23        | 9630TA26        | 0690103D003CER000QZ |
| Edwards Certitude introducer sheath set          | 9620IS18        | 9620IS21        |                 | 0690103D003CIS000SL |
| inflation device                                 | 96402           |                 |                 | 0690103D003IND000TG |
| Edwards crimper                                  | 9600CR          |                 |                 | 0690103D003CRI000TH |

## 14.0 Expected Lifetime of the Device

The Edwards transcatheter heart valve has been subjected to rigorous pre-clinical durability testing per the valve testing requirements and in clinical studies and post market studies. The valves were successfully tested to 5 years of simulated wear. In addition, clinical data show

durability with follow-up to 5 years. The actual lifetime performance is continuing to be studied and varies from patient to patient.

## 15.0 Patient Information

A patient implant card is provided with each THV. After implantation, please complete all requested information and provide the implant card to the patient. The serial number is found on the package. This implant card allows patients to inform healthcare providers what type of implant they have when they seek care.

## 16.0 Recovered THV and Device Disposal

The explanted THV should be placed into a suitable histological fixative such as 10% formalin or 2% glutaraldehyde and returned to the company. Refrigeration is not necessary under these circumstances. Contact Edwards Lifesciences to request an Explant Kit.

Used devices may be handled and disposed of in the same manner as hospital waste and biohazardous materials. There are no special risks related to the disposal of these devices.

## 17.0 Clinical Studies

Refer to the SSCP for clinical benefits.

## 18.0 References

1. Bapat V, Attia R, Thomas M. Effect of Valve Design on the Stent Internal Diameter of a Bioprosthetic Valve: A Concept of True Internal Diameter and Its Implications for the Valve-in-Valve Procedure. JACC: Cardiovascular Interventions. Vol. 7, No. 2. 2014: 115-127.

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
|                                                                                   |                                                                                             |                                                                                             |                                                                                             |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------|
|  | Non-pyrogenic                                                                               | Ikke-pyrogen                                                                                | Pyrogeeniton                                                                                |
| <b>MD</b>                                                                         | Medical device                                                                              | Medisinsk utstyr                                                                            | Lääkinnällinen laite                                                                        |
|  | Contains biological material of animal origin                                               | Inneholder biologiske materialer av animalsk opprinnelse                                    | Sisältää eläinperäistä biologista materiaalia                                               |
| <b>Rx only</b>                                                                    | Caution: Federal (USA) law restricts this device to sale by or on the order of a physician. | Caution: Federal (USA) law restricts this device to sale by or on the order of a physician. | Caution: Federal (USA) law restricts this device to sale by or on the order of a physician. |
|  | Time & Temperature Sensitive                                                                | Tids- og temperaturfølsom                                                                   | Aika- ja lämpöherkkä                                                                        |
|  | Contains hazardous substances                                                               | Inneholder farlige stoffer                                                                  | Sisältää vaarallisia aineita                                                                |
| <b>SZ</b>                                                                         | Size                                                                                        | Størrelse                                                                                   | Koko                                                                                        |

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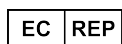2023-06  
10055170001 C  
© Copyright 2023, Edwards Lifesciences LLC  
All rights reserved.

**Edwards Lifesciences LLC**  
One Edwards Way  
Irvine, CA 92614 USA

**Edwards Lifesciences B.V.**  
Verlengde Poolseweg 16  
4818 CL Breda, Netherlands

Telephone +1 949.250.2500  
+1 800.424.3278  
FAX +1 949.250.2525

Web IFU

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**Edwards Lifesciences GmbH**